

FINALLY!

*A portable device you can configure
for all your monitoring needs*

petMAP+ II

Multi-Parameter Monitor



Developed by:
**RAMSEY
MEDICAL INC**

Manufactured,
Distributed and
Serviced by:



www.cardiocommand.com

**Portability, Accuracy, Flexibility
and Simplicity like no other!**

www.petmap.com



petMAP+ II

Multi-Parameter Monitor

Exam Room & OR Monitoring With One Handheld Device

Buy a base unit . . .

7400	NIBP/Temp Device
7402	NIBP/SpO2/Temp Device
7404	NIBP/ <i>NONIN</i> SpO2/Temp Device

petMAP+II provides:

- Proven petMAP BP technology with PPO (PetMAP Proprietary Optimizations). Species and cuff site settings for improved accuracy
- O2 saturation option with pleth display
- Temperature monitor using rectal or esophageal probe
- Plugin CO2 and ECG modules available
- Easy to use "no menu" touch screen display
- Real-time display of BP oscillometric "envelope"
- Selectable automatic BP cycle times
- SD card option stores trend data & screen shots
- NiMH batteries and wall charger
- Optional pole mount and cover



Blood Pressure Range: Systolic BP: 40 - 265 mmHg
MAP: 30 - 200 mmHg
Diastolic BP: 25 - 195 mmHg

Heart Rate Range: 40 - 240 bpm (both BP and SpO2)

SpO2: +/- 3% in 80 - 100% range

Temperature: +/- .2° C in 29° C - 46° C range

Dimensions: 6.0"H x 3.2"W x 1.25"D

Weight: 1 lb (with NiMH batteries)

Battery : Std: four AA NiMH, Optional: Alkaline


Then add optional modules.

Just Plug and Play!

petMAP+ II

ECG
Module

- Displays ECG waveform and HR
 - 3 electrode system; lead II for HR
 - Esophageal catheters available.
- 2 sizes. Displays ECG, HR and temperature

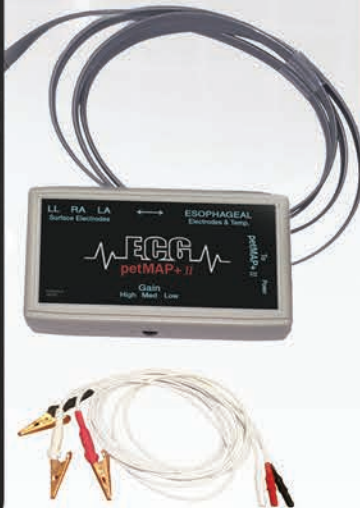

REF#: 9026

Heart Rate Range: 20-260 BPM

Heart Rate Accuracy: 3% +1 digit

petMAP+ II

ETCO₂ Module

- Mainstream system
- Provides capnograph, CO2 numerics (ETCO2, Insp CO2, Rt CO2) and Resp Rate from CO2.
- Includes large and small airway adapters



REF#: 9027

Heart Rate Range: 20-260 BPM

Heart Rate Accuracy: 0-40 ± 2 mmHg
41-70 ± 5 % of reading
71-100 ± 8 % of reading

RR Accuracy: ± 2%; ± 1 BPM

RR Range: 40-140 BPM

petMAP+ II

Multi-Parameter Monitor

CritterCuffs™



CritterCuffs included with petMAP+ II System

Blood Pressure Cuffs

Seven sizes included with petMAP+ II

Size	Item no.
2.0cm	8050
2.5cm	8051
3.0cm	8055
3.5cm	8052
4.0cm	8056
4.5cm	8053
5.5cm	8054
5.0cm	8067*
6.5cm	8057*
8.0cm	8058*
10.0cm	8059*
13.0cm	8060*



*Sold separately, not included in petMAP+ II System

RamseyCuffs™



*Sold separately, not included in petMAP+ II System

Available Accessories

Pole Mount

Item #: 9029



Cover

Item #: 9028 (Black, Grey, Blue or Red)



Rectal Temperature Probe

Item #: 8041



Esophageal Temperature Probe

Item #: 8040



Temperature Extension Cable

Item #: 9014



International Adaptor Plugs

Item #: 9008



Cat Catheter (ECG & Temp.)

Item #: 3205-0001



Requires ECG Module

Dog Catheter (ECG & Temp.)

Item #: 3007-0001



Requires ECG Module

3 Lead Wire Set

Item #: 9030



Optional Nonin SpO2 Sensors



Item #: 9018
Transflectance Sensor

Item #: 9019
Small Flex Sensor

CO2 Large Airway Adapter

Item #: 9034



CO2 Small Airway Adapter

Item #: 9035



CONTACT US TODAY TO PLACE YOUR ORDER

(800) 231-6370

(813) 289-5555

Take advantage of our 30 day satisfaction guarantee.

For more information visit www.petmap.com

Manufactured,
Distributed and
Serviced by:



Developed by: **RAMSEY MEDICAL INC**
R400304-01
09/15/16